



Surgical perfection. For life.

Setting new low-dose standards in the OR

The STILLE imagiQ2™ surgical imaging table represents an important advance in developing the low-dose technology field, making the OR an even safer place, for patients and for the entire surgical team.


More efficient and shorter procedure times

The True Free Float® technology* of the STILLE imagiQ2 gives unprecedented control and freedom of movement in any desired direction while providing state-of-the-art smoothness like an extension of the surgeon's hand, allowing for shorter procedure times and higher efficiency in the OR.

Choose a low-dose enabler

The unique high-performance carbon fiber table top*, provides unbeatable transparency of 0.4 mm Al equivalency, enabling the C-arm to reduce the radiation dose without compromising image quality.

*Patent pending

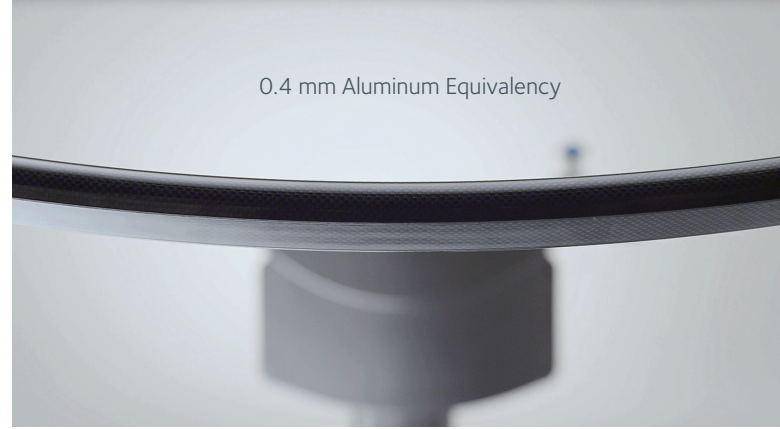

Optimizing C-arm capabilities

With its ultra-slim design and super-low base, the imagiQ2 enhances C-arm access, making it the ideal C-arm complement for the mobile Hybrid OR or office suite. Furthermore, the low dose table top of the Stille imagiQ2 contributes to avoiding C-arm overheating.

Continuous image focus


The patented technology of the iso-roll cradling motion ($\pm 15^\circ$) helps maintain image focus even while in full lateral roll, minimizing patient and C-arm repositioning, as well as reducing scattered radiation.

STILLE imagiQ2

- Translucent length: patients as tall as 6'5" / 196 cm can be fully imaged with headrest (5'6" / 169.2 cm without headrest). 360° unobstructed translucent length: 3'6" / 107.2 cm with headrest (2'7" / 80.4 cm without headrest).
- Patient weight capacity: 225 Kg (496 lbs). Can be upgraded to handle patients as heavy as 300 Kg (661 lbs) with the *Load Upgrade Pack*.
- True Free Float® allows smooth longitudinal (35.4" / 90 cm), lateral (9.9" / 25 cm) and diagonal (36.8" / 93.5 cm) movements.
- Upgrades to fit your needs: Load upgrade pack, performance pack, mobility enhancement pack.

Learn more at: www.stille.se



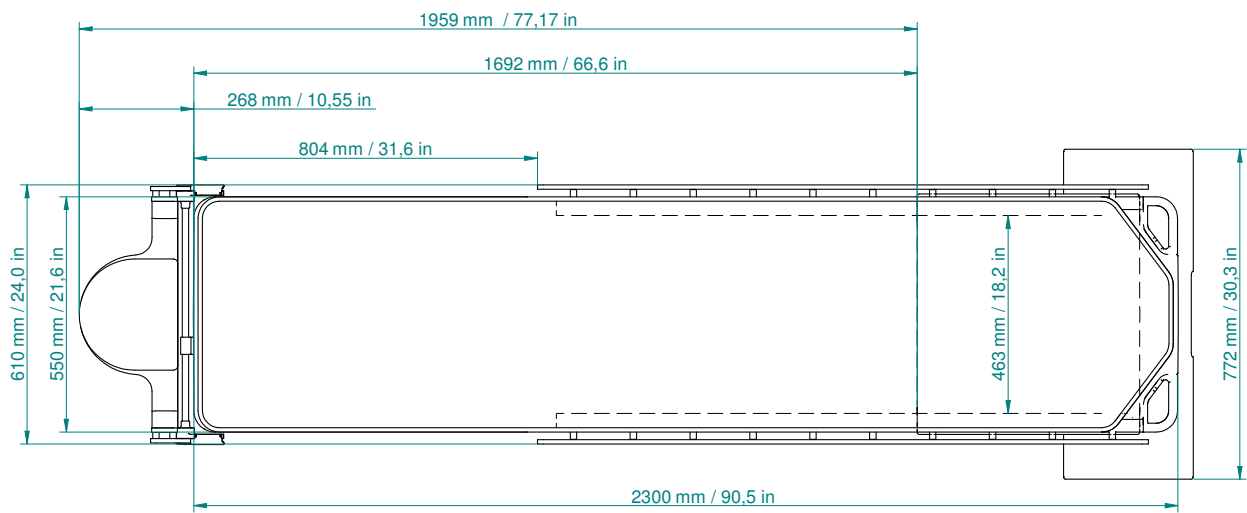

Product specifications of the Stille imagiQ2

Feature	Description
Patient weight	496 lb / 225 kg and 661 lb / 300 kg (with load upgrade pack)
Table top length , width	L= 90.5" / 230 cm ; W= 21.6" / 55 cm
Translucent length	77.17" / 196 cm with headrest (66.6" / 169.2 cm without headrest)
Translucent table top	Proprietary, patent pending carbon fiber structure. Al 0.4
Float lock	Stille quick-lock system
Float (longitudinal, lateral, diagonal)	35.4" x 9.9" x 36.8" / 90 cm x 25 cm x 93.4 cm
Operation conditions	100% continuous
AC power: 100-250 Volts	Direct AC power + battery DC power
Length, width, height	L = 92" / 234 cm; W = 30.3" / 77.2 cm; H = 28"- 43" / 71-109 cm
Trendelenburg	±25°
Lateral roll (iso-roll)	±15°
Wheels and steering	4 swiveling back wheels
Side rails	Fixed along 2/3 of the table top
Exchangeable table top	Yes, future-proof platform for extended use
Transport handle	Yes (Color: Black for 225-kg version, silver for 300-kg version)

Customizing the imagiQ2 to fit your needs

Upgrade features	Description
Mobility enhancement pack	When the table is used for frequent intra-clinic transportation, use the mobility enhancement pack for improved mobility. All wheels are designed to swivel and an additional fifth wheel is included to provide optimal mobility. Ref: 535-1711
Load upgrade pack	Maximize the capabilities of the imagiQ2 with increased patient weight capacity up to 661 lb / 300 kg. Ref: 535-1710
Performance pack	For maximum uptime, choose the performance pack. It includes a high-performance long-term battery pack with up to 1 week of durability. Ref: 535-1712