



Surgical perfection. For life.

Setting new low-dose standards in the OR The imagiQ2 is the Low-Dose Enabler™

The STILLE imagiQ2™ surgical imaging table represents an important advance in developing the low-dose technology field, making the OR an even safer place, for patients and for the entire surgical team.

More efficient and shorter procedure times

The True Free Float® technology* of the STILLE imagiQ2 gives unprecedented control and freedom of movement in any desired direction while providing state-of-the-art smoothness like an extension of the surgeon's hand, allowing for shorter procedure times and higher efficiency in the OR.

Choose a Low-Dose Enabler™

The exceptional high-performance carbon fiber table top*, provides outstanding transparency of 0.4mm Al equivalency, enabling the C-arm to reduce the radiation dose without compromising image quality. That is 60% more translucency than average tables on the market.

*Patent pending



MED TECH TOP 10
OUTLOOK
MEDICAL IMAGING
SOLUTION PROVIDERS - 2018

Optimizing C-arm capabilities

With its ultra-slim design and super-low base, the imagiQ2 enhances C-arm access, making it the ideal C-arm complement for the mobile Hybrid OR or office suite. These features minimize the need for C-arm adjustments that could contribute to overheating.

Continuous image focus

The patented technology of the Iso Roll™ cradling motion ($\pm 15^\circ$) helps maintain image focus even while in full lateral roll, minimizing patient and C-arm repositioning, as well as reducing scattered radiation.

STILLE imagiQ2

- Translucent length: patients as tall as 6'5" / 196 cm can be fully imaged with headrest (5'6" / 169.2 cm without headrest). 360° unobstructed translucent length: 3'6" / 107.2 cm with headrest (2'7" / 80.4 cm without headrest).
- Patient weight capacity: 661 lbs. (300 kg).
- True Free Float® allows smooth longitudinal (35.4" / 90 cm), lateral (9.9" / 25 cm) and diagonal (36.8" / 93.5 cm) movements.

Gold-standard service and warranty STILLE is dedicated to providing gold-standard customer service. We believe in our products: imagiQ2 is backed with a 2-year warranty on labor and 3 years on parts.

Now Featuring

Standard Weight Capacity **661lbs/300kg**

Extended battery life- Up to five days*

*battery life contingent upon pt. weight and table use

Learn more at: www.stille.se

Product specifications of the STILLE imagiQ2

| Feature | Description |
|---|---|
| Patient weight | 661 lb/300 kg |
| Table top length , width | L= 90.5" / 230 cm ; W= 21.6" / 55 cm |
| Translucent length | 77.17" / 196 cm with headrest (66.6" / 169.2 cm without headrest) |
| Float (longitudinal, lateral, diagonal) | 35.4" x 9.9 " x 36.8" / 90 cm x 25 cm x 93.4 cm |
| AC power: 100-250 Volts | Direct AC power + Power Performance Battery DC power |
| Length, width, height | L = 92" / 234 cm; W = 30.3" / 77.2 cm; H = 28"- 43" / 71-109 cm |
| Trendelenburg | ±20° |
| Lateral roll (Iso Roll™) | ±15° (Iso Roll™ keeps image focus & reduces scattered radiation) |
| Side rails | Fixed along 2/3 of the table top |

Standard configuration includes: 2 hand controls, 1 extension cable, 1 pan handle, mattress.

Available Accessories

(US)


| Accessory | Part number |
|----------------------------------|-------------|
| Long detachable side rail | 535-1720 |
| Short detachable side rail | 535-1724 |
| Fistula arm board | 535-1730 |
| Head rest, adjustable | 535-1742 |
| Catheter tray 1000 mm | 535-1750 |
| Catheter tray 1500 mm | 535-1751 |
| Catheter tray 500 mm | 535-1752 |
| Foot Control | 535-1762 |
| Basic armboard | 512-5000 |
| Anesthesia screen | 512-5030 |
| IV pole w/o clamp | 512-5010 |
| IV pole clamp | 512-5013 |
| Foot board carbon fiber with pad | 512-5040 |
| Deluxe patient strap | O-RSX2 |
| Radiation shield | 535-1773 |



Radiation Shield



Basic arm board

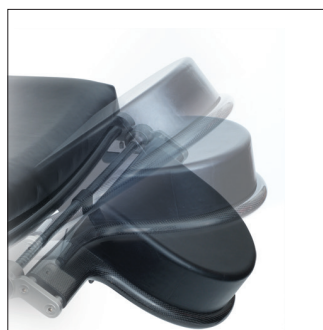

Adjustable procedure headrest



Detachable short side rail

Accessory Package Available Part 535-1790

| Accessory package includes: | Part number |
|-----------------------------|-------------|
| Fistula arm board | 535-1730 |
| Catheter Tray 39" | 535-1750 |
| Radiation Shield (grey) | 535-1772 |
| Long Side Rail | 535-1720 |
| Short Side Rail | 535-1724 |
| Deluxe Patient Strap | O-RSX2 |

Package offers 20% savings over individual items purchased separately.



Translucent fistula arm board