



Surgical perfection. For life.



Best Surgical Imaging Table

Shortening Procedure Times while Reducing Radiation Exposure

The STILLE imagiQ2™ surgical imaging table represents an important advance in developing the low-dose technology field, making the OR an even safer place for patients and for the entire surgical team.

Increase Efficiency and Shorten Procedure Times

The True Free Float™ technology of the STILLE imagiQ2 gives unprecedented control and freedom of movement in any desired direction while providing state-of-the-art smoothness. Like an extension of the surgeon's hand, the True Free Float™ increases efficiency in the OR and shortens procedure times.

The patented technology of the ISO ROLL™ cradling motion ($\pm 15^\circ$) helps maintain image focus even while in full lateral roll, minimizing patient and C-arm repositioning, as well as reducing scatter radiation.




Reduce Radiation Exposure and Enhance Image Quality

The unique high-performance carbon fiber table top provides outstanding radiolucency of 0.4 mm Al equivalency, enabling the C-arm to reduce the radiation exposure without compromising image quality. That is up to 60% more radiolucent than average tables on the market.

Its ultra-slim design, super low-base and retractable wheels, maximize C-arm access and allow the image receiver to come closer to the patient, thus reducing radiation exposure while enhancing image quality.


Choose a Low-Dose Enabler™

Taken together, the True Free Float™ technology, outstanding table top radiolucency and retractable wheels, is what makes the STILLE imagiQ2™ a true Low-Dose Enabler™ and the ideal C-arm complement for the mobile Hybrid OR or office suite.

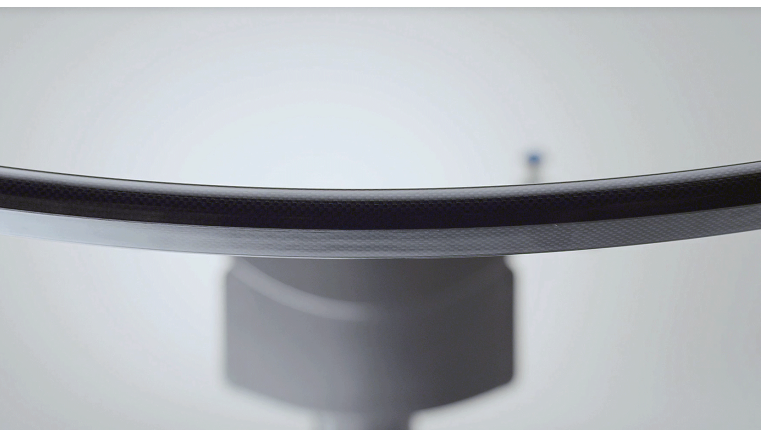

Low-dose table top with rounded edges and the highest radiolucency in the world for floating tables. Specifically designed for oblique imaging and eliminating image artifacts while enhancing image quality.

Learn more at: www.stille.se/imagiq2

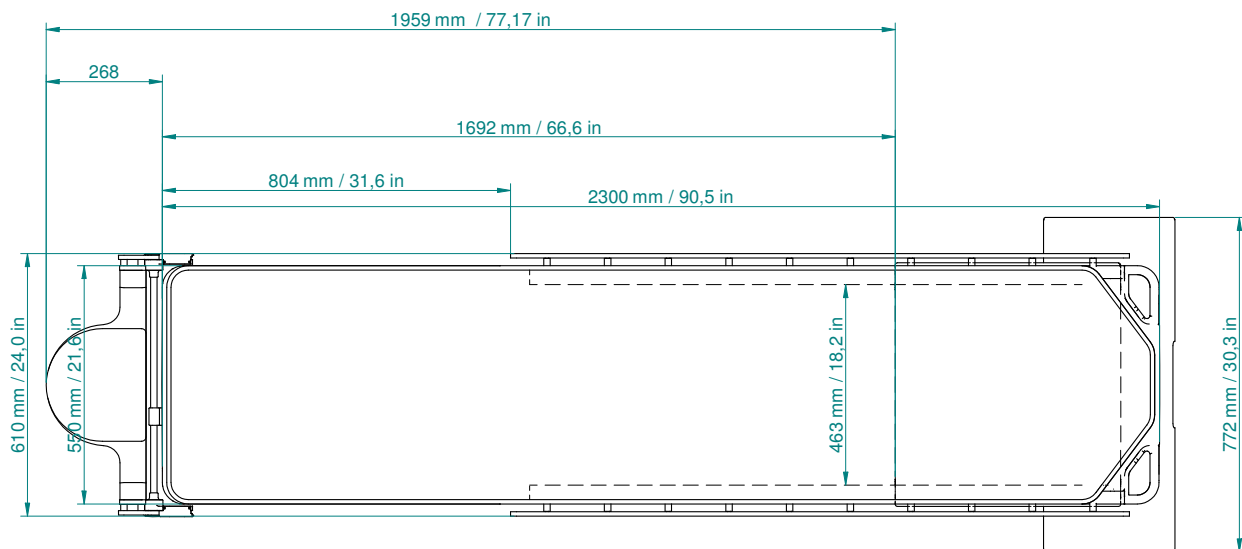

Product Specifications of the STILLE imagiQ2™

Feature	Description
Weight capacity	496 lb / 225 kg (optional 661 lb / 300 kg with load upgrade pack)
Table top length , width	L= 90.5" / 230 cm ; W= 21.6" / 55 cm
Radiolucent length	77.17" / 196 cm with headrest (66.6" / 169.2 cm without headrest)
Radiolucent table top	Proprietary, carbon fiber structure. Al 0.4 mm
Float & lock	STILLE True Free Float® and STILLE quick-lock system
Float (longitudinal, lateral, diagonal)	35.4" x 9.9" x 36.8" / 90 cm x 25 cm x 93.4 cm
Operation conditions	100% continuous
AC power: 100-250 Volts	Direct AC power + battery DC power
Length, width, height	L = 92" / 234 cm; W = 30.3" / 77.2 cm; H = 28"- 43" / 71-109 cm
Trendelenburg	±25° (with load upgrade pack ±20°)
Lateral roll (ISO ROLL™)	±15°
Wheels and steering	4 swiveling back wheels
Side rails	Fixed along 2/3 of the table top
Detachable table top	Yes, future-proof platform for extended use and easy transport.
Transport handle	Yes (Color: dark grey for 225 kg-version, silver for 300 kg-version)

Customizing the imagiQ2™ to Fit Your Needs

Upgrade Feature	Description
Mobility Enhancement Pack	When the table is used for frequent intra-clinic transportation, use the mobility enhancement pack for improved mobility. All wheels are designed to swivel and an additional fifth wheel is included to provide optimal mobility. Ref: 535-1711
Weight Upgrade Kit	Maximize the capabilities of the imagiQ2™ with increased weight capacity up to 661 lb / 300 kg. Ref: 535-1710



Best Surgical Imaging Table



Stille AB Ekbacken 11, SE-644 30 Torshälla, Sweden
Phone: +46 8 588 580 00 • sales@stille.se • www.stille.se