




STILLE Medstone

Surgical Imaging Tables

Superior versatility

for a wide range of
interventional procedures

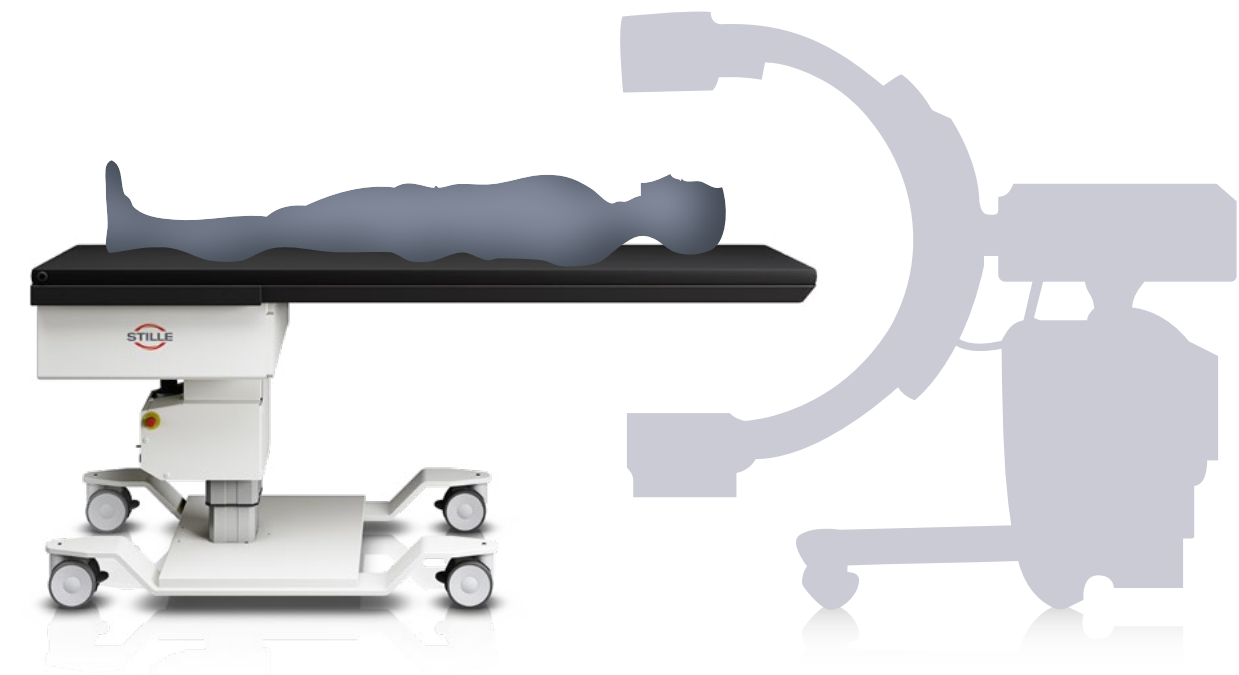

SUPINE
POSITION



PRONE
POSITION



LATERAL
POSITION



SPINAL
POSITION



TRENDELENBURG
POSITION



REVERSE
TRENDELENBURG
POSITION



- ✓ Affordable quality
- ✓ Intuitive and user-friendly
- ✓ Large imaging areas

The vast selection of Medstone models support an extensive array of surgical interventions. The quality and reliability of our C-arm surgical imaging tables contribute to enhancing efficiency during procedures at an economical price.

Wide range of movement for precise patient positioning. Whether motorized or via the free floating function – optimal patient positioning is ensured intuitively and quickly.

Up to five motorized movements available



Height adjustment

For easy patient transition and mobility as well as adjustment to the physician's height. With an elevation of up to 108 cm, the table allows the patient to be positioned close to the image receptor which increases image quality.



Longitudinal and lateral travel

Motorized movement enables precise longitudinal and lateral positioning for X-ray imaging as the C-arm does not have to be moved to capture anatomy of interest.



Trendelenburg

In situations where the patient's head needs to be raised or lowered.



Lateral tilt

Minimizes the need for repositioning the C-arm during procedures.



Free Float technology with Medstone3 CV

Allows complete freedom of movement horizontally and diagonally for more reliability and control when operating the table. This enables excellent accuracy and the reduction of procedure time. Ideal for applications such as following contrast media during endovascular procedures (bolus chase).

Over 1,000 units sold worldwide

STILLE Medstone tables are reliable in multiple applications. Customers appreciate the intuitive and user-friendly operation. In combination with the carbon fiber tabletop, the Medstone tables contribute to high image quality.



STILLE Medstone PM

- ▶ offers up to five motorized movements
- ▶ available with different table tops: contoured with a facial cutout for prone patient positioning and rectangular without a facial cutout for supine patient positioning
- ▶ for every demand: Medstone2 PM is the entry-level model, while Medstone5 PM is ideal for complex treatments
- ▶ 227 kg weight capacity

Standard configuration
mattress, three patient straps, hand control



- Contoured:
- ✓ Orthopedics
 - ✓ Pain management
- Rectangular:
- ✓ Electrophysiology
 - ✓ Peripheral vascular
 - ✓ Veterinary
 - ✓ Multi-disciplinary

	Medstone2 PM	Medstone3 PM	Medstone3 PM Rectangular	Medstone4 PM	Medstone5 PM
Part number	500-1011	500-1005	500-1043	500-1004	500-1001
Table top type	contoured ⁽¹⁾	contoured ⁽¹⁾	rectangular	contoured ⁽¹⁾	contoured ⁽¹⁾
Table top measurements (cm)	56 x 208	56 x 208	56 x 208	56 x 208	56 x 208
360° imaging length (cm)	147	147	147	147	147
Fixed side rail (cm)	no	no	no	no	30.5
Infectious control covers	no	no	no	no	yes
Height adjustment ⁽²⁾ (cm)	✓ 65.3 – 100.6	✓ 62.6 – 97.6	✓ 62.6 – 97.6	✓ 73.0 – 108.0	✓ 65.5 – 101.0
Trendelenburg ⁽³⁾ (± 15°)	✗	✓	✓	✓	✓
Lateral tilt ⁽³⁾ (± 15°)	✓	✗	✗	✓	✓
Longitudinal travel ⁽²⁾ (40 cm)	✗	✓	✓	✓	✓
Lateral travel ⁽²⁾ (15 cm)	✗	✗	✗	✗	✓

⁽¹⁾ with facial cutout

⁽²⁾ ± 1.25 cm tolerance

⁽³⁾ ± 2.0° tolerance

STILLE Medstone ERCP

- ▶ offers up to five motorized movements
- ▶ ideal for procedures with the patients in lateral position and/or for larger patients
- ▶ comes with an extra wide table top to easily facilitate abdominal fluoroscopy and endoscopic procedures
- ▶ 227 kg weight capacity

Standard configuration
mattress, three patient straps, hand control



- ✓ Abdominal
- ✓ Gastroenterology
- ✓ Bronchoscopy
- ✓ Orthopedics
- ✓ Fluoroscopic procedures

	Medstone2 ERCP	Medstone3 ERCP	Medstone3 ERCP	Medstone5 ERCP
Part number	500-1041	500-1007	500-1044	500-1031
Table top type	wide	wide	wide	wide
Table top measurements (cm)	66 x 208	66 x 208	66 x 208	66 x 208
360° imaging length (cm)	147	147	147	147
Fixed side rail (cm)	no	no	no	no
Infectious control covers	no	no	no	yes
Height adjustment ⁽¹⁾ (cm)	✓ 66.5 – 101.9	✓ 63.5 – 98.8	✓ 63.9 – 98.9	✓ 67.0 – 102.0
Trendelenburg ⁽²⁾ (± 15°)	✗	✗	✓	✓
Lateral tilt ⁽²⁾ (± 15°)	✓	✗	✗	✓
Longitudinal travel ⁽¹⁾ (40 cm)	✗	✓	✓	✓
Lateral travel ⁽¹⁾ (15 cm)	✗	✓	✗	✓

⁽¹⁾ ± 1.25 cm tolerance

⁽²⁾ ± 2.0° tolerance

STILLE Medstone CV

- ▶ the STILLE Medstone3 CV table with free floating table top for head-to-toe imaging with complete accuracy, speed and freedom of movement
- ▶ the STILLE Medstone3 CV has a float control unit with magnetic locking release as well as foot and hand controls
- ▶ with the STILLE Medstone5 CV table, the series offers five motorized movements, including Trendelenburg
- ▶ 227 kg weight capacity

Standard configuration
mattress, three patient straps, hand control;
Medstone3 CV also includes a foot control and a pan handle



- ✓ Cardiovascular
- ✓ Interventional Radiology
- ✓ Endovascular
- ✓ Fluoroscopic procedures
- ✓ Vascular
- ✓ Electrophysiology (Pacemaker, ...)

	Medstone3 CV	Medstone5 CV
	with free float	
Part number	500-1025	500-1002
Table top type	rectangular	rectangular
Table top measurements (cm)	61 x 208	56 x 208
360° imaging length (cm)	76 ⁽¹⁾	147
Fixed side rail (cm)	93.0	30.5
Infectious control covers	no	yes
Height adjustment ⁽²⁾ (cm)	✓ 67.6 – 107.6	✓ 65.5 – 101.0
Trendelenburg ⁽³⁾ (± 15 °)	✗	✓
Lateral tilt ⁽³⁾ (± 15 °)	✗	✓
Longitudinal travel ⁽²⁾ (cm)	✓ 91.0	✓ 40.0
Lateral travel ⁽²⁾ (cm)	✓ 20.0	✓ 15.0

⁽¹⁾ total imaging length: 167 cm

⁽²⁾ ± 1.25 cm tolerance

⁽³⁾ ± 2.0° tolerance

Recommended accessories. Here is a selection out of the wide range of accessories available. The following is just a selection of the expansive offerings.

Accessory	Type	Part number	Compatible with
X-ray protective kit	grey blue	535-1770 535-1771	all models all models
X-ray protective loop	grey blue	535-1772 535-1773	all models all models
Foot control	5-movement 4-movement 3-movement 3-movement	505-1109 505-1108 505-1139 505-1138	Medstone5 PM / CV / ERCP Medstone4 PM Medstone3 ERCP (500-1007 only) Medstone3 PM / ERCP (500-1005 and -1044 only)
Detachable side rails	2 x 80.5 cm, symmetrical 2 x 80.5 cm, symmetrical 2 x 112.5 cm, non-symmetrical 2 x 112.5 cm, non-symmetrical	505-1179 505-1181 505-1175 505-1177	Medstone2, 3 and 4 PM (500-1004 only), Medstone3 CV Medstone2 and 3 ERCP (500-1040 and -1044 only) Medstone PM, Medstone5 CV Medstone ERCP
Arm board	snap-on arm board radiolucent fistula arm board adjustable arm board	505-1111 505-1118 514-50-3SR025	all models all models excl. Medstone ERCP all models
Side gates	foldable, pair	514-100530	Medstone ERCP
Catheter tray	68.3 cm	505-1131	Medstone3 CV
Accessory package	vascular, incl. 505-1179, 535-1772, 505-1118 ERCP incl. 505-1181, 514-50-3SR025, 514-100530 ERCP incl. 505-1181, 514-100530	505-1188 505-1189 505-1190	Medstone3 CV Medstone ERCP Medstone ERCP

Please see the accessory brochure for the full list.

STILLE – Surgical perfection. For life.

STILLE is one of the oldest medical device companies in the world. At STILLE, we develop, produce and sell high-quality surgical instruments and surgical imaging tables to physicians worldwide, and have done so for about 180 years.

STILLE's offering consists of a strong brand and products renowned for their superior quality. We manage our business in two distinct business areas – Surgical Instruments and Surgical Imaging Tables.

Surgical Instruments

In our Surgical Instruments business area, we have a complete offering of handcrafted, forged instrumentation for all fields of open surgery, with a special emphasis on cardiovascular surgery, and reconstructive and aesthetic plastic surgery.

Our instruments feature delicate design, industry leading durability and a unique feel, which has made STILLE a well-known and highly respected supplier by surgeons throughout the world.

Surgical Imaging Tables

In our Surgical Imaging Tables business unit, we offer specialty tables for minimally invasive procedures, with pricing, safety and features optimized for both in- and outpatient settings. The STILLE Medstone mobile surgical imaging tables contribute to enhance efficiency in the OR, all while representing valued priced solutions with both quality and reliability.

Learn more at www.stille.se

STILLE AB

Ekbacken 11

644 30 Torshälla

Sweden

☎ +46 8 588 580 00

✉ sales@stille.se

🌐 www.stille.se



Surgical perfection. For life.

These products may not be available in all regions. Please contact your STILLE-representative for more information on market availability.

This document correctly describe these devices at the time the document was printed. Please always refer to the user manual for any product amendments.

The product images shown are for illustration purposes only and may not be an exact representation of the product.