




Surgical perfection. For life.

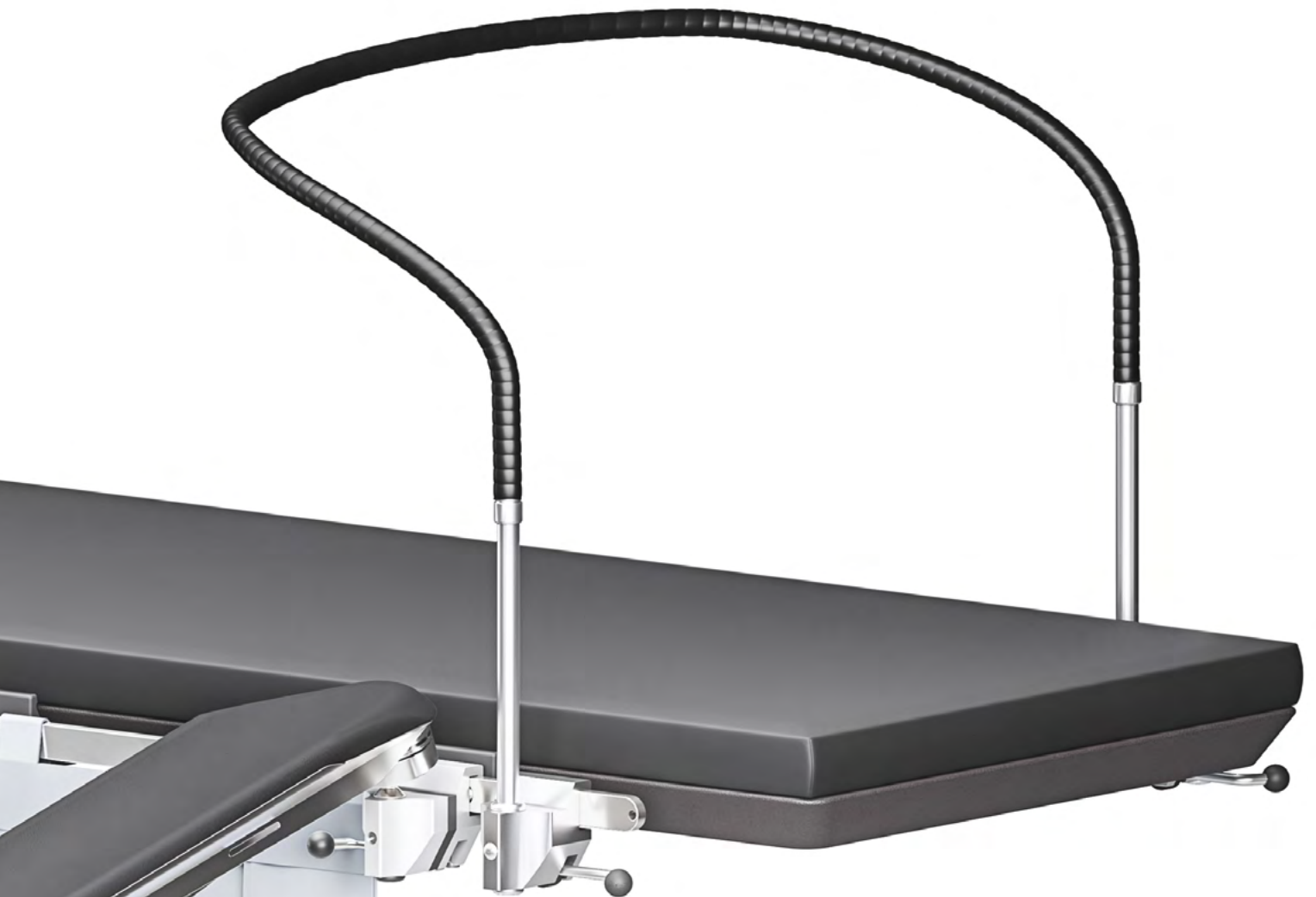

# STILLE Medstone

## Accessory catalog



# Table of contents

<b>Introduction</b>	<b>4</b>
<b>Accessory packages</b>	<b>6</b>
<b>Detachable side rails</b>	<b>8</b>
70.0 cm / 27.5" detachable side rails	8
80.5 cm / 31.7" detachable side rails	9
112.5 cm / 44.3" detachable side rails	10
Accessory table clamps	11
<b>Anesthesia accessories</b>	<b>12</b>
Anesthesia screen	12
IV pole	12
<b>Arm, leg &amp; foot supports</b>	<b>14</b>
Arm boards	14
Arm & leg guards	16
Tibia support	17
<b>Control devices</b>	<b>18</b>
Extra pan handle	18
Extension cable	18
Foot controls	19
<b>General body supports &amp; barriers</b>	<b>22</b>
Lateral support	22
Side gates	22
Patient restraints / Leg & body strap	23
<b>Radiation protection</b>	<b>24</b>
<b>Surgical table extensions</b>	<b>26</b>
Catheter tray	26
Instrument tray	26
Widening sections	27
<b>Index</b>	<b>28</b>

STILLE Medstone Surgical table accessories offer an extensive range of solutions providing support, patient safety and flexible positioning to OR procedures. All of our accessories are carefully designed and made with strong emphasis on quality.

How to navigate through our catalogue?  
All stored into one document for easy reference, STILLE Medstone accessories have been grouped and presented individually as illustrated below.

Page | CATEGORY OF THE ITEM

CATEGORY OF THE ITEM

Product name or specific characteristic(s)

Description of the item: use, characteristics, quick mounting instructions, material, etc.

✔ **Included:**  
Parts and numbers related to the item

⚠ **Required:**  
Other accessories needed to mount or to use the item

⊕ **Optional:**  
Alternative or suggestion of combination of items for a smoother workflow or optimal use

Compatibility

✔ Medstone CV  
✔ Medstone ERCP  
✔ Medstone PM

Order information

EU: XX-XX-XXXX  
US: XX-XX-XXXX


Illustration of the item. May not be an exact representation of the product.

List of compatible tables with the item. It can be all models of a type of Medstone (CV, PM or ERCP) or, limited to few references only (e.g. Medstone3 PM 500-1005)

Part number\* of the item for your order. Part number might be different depending on the region. Most regions outside the USA use the EU version. Products may not be available in all regions. Please contact your STILLE-representative for more information on market availability.

\* To search by part number, simply refer to the index at the end of the catalogue.

ACCESSORY PACKAGES

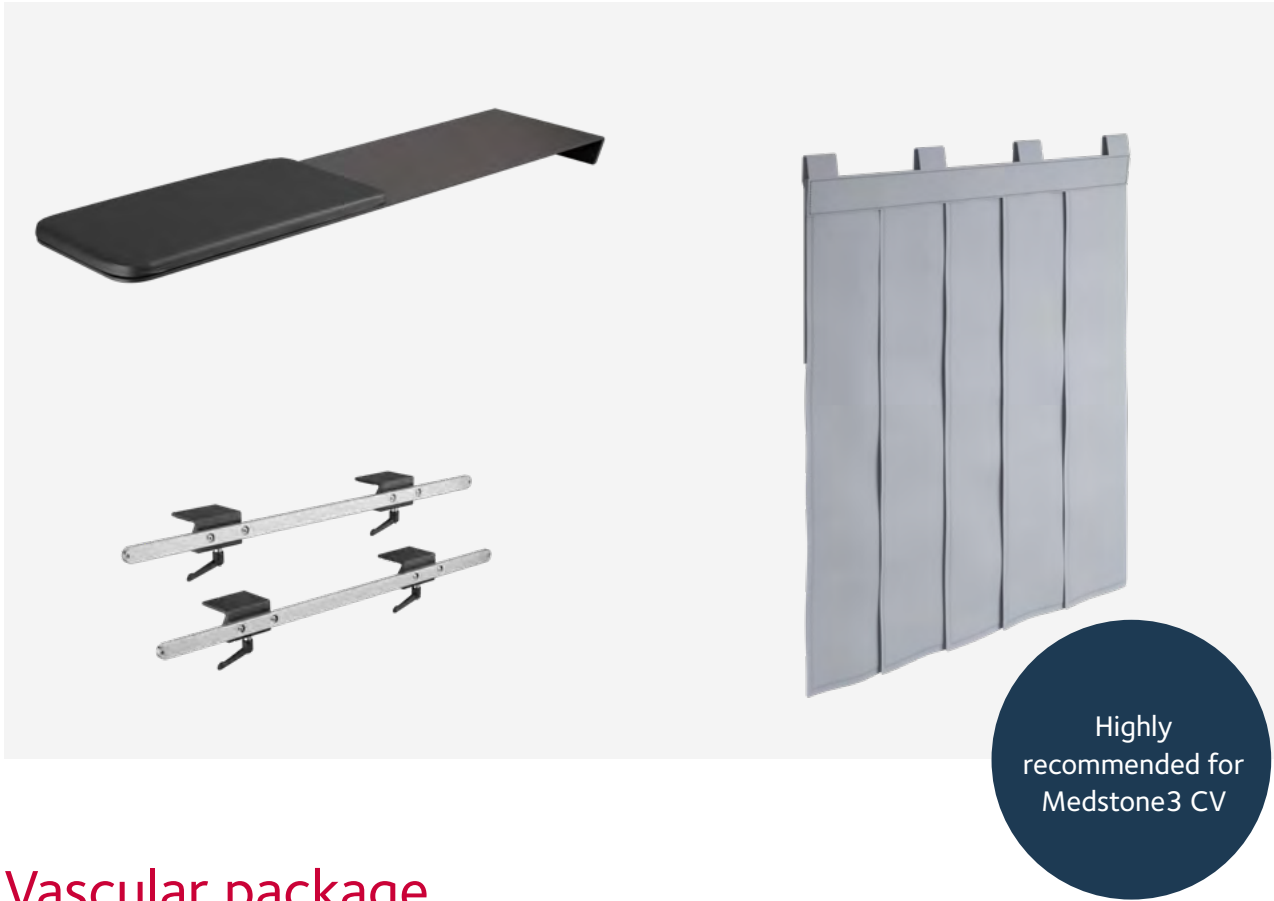

Vascular package  
Medstone3 CV

A procedure pack consisting of side rails, arm board, radiation shield, all essential accessories for using Medstone3 CV. Not available in the USA as a package; items can be ordered individually.

- ✓ **Included:**  
Pair of 80.5 cm detachable side rails  
Fistula arm board with a mattress  
Radiation shield loop grey
- ! **Required:**  
-
- + **Optional:**  
-

Compatibility

- ✓ Medstone3 CV
- ✓ Medstone2 PM
- ✓ Medstone3 PM
- ✓ Medstone4 PM (500-1004 only)

Order information

**EU:** 505-1188  
**US:** Not available

Included accessories

<b>Fistula arm board</b> page 14 <b>EU:</b> 505-1118 <b>US:</b> 505-1118	<b>Radiation shield loop grey</b> page 24 <b>EU:</b> 535-1772 <b>US:</b> 535-1772
<b>Pair of side rails</b> <b>80.5 cm / 31.7" length</b> page 09 <b>EU:</b> 505-1179 <b>US:</b> 505-1180	



ERCP package  
Medstone

A procedure pack consisting of side gates to secure patients, side rails, and a arm board. All essential accessories for ERCP procedures. Not available in the USA as a package; items can be ordered individually.

- ✓ **Included:**  
Pair of 80.5 cm detachable side rails  
Pair of side gates  
Adjustable arm board
- ! **Required:**  
-
- + **Optional:**  
ERCP package without the adjustable armboard:  
**EU:** 505-1190  
**US:** Not available

Compatibility

- ✓ Medstone ERCP

Order information

**EU:** 505-1189  
**US:** Not available

Included accessories

<b>Adjustable armboard</b> <b>compl. DIN</b> page 14 <b>EU:</b> 514-50-3SR025 <b>US:</b> 514-50-3SR027	<b>Side gates, foldable / pair</b> page 22 <b>EU:</b> 514-100530 <b>US:</b> 514-100530-US
<b>Pair of side rails (W)</b> <b>80.5 cm / 31.7" length</b> page 09 <b>EU:</b> 505-1181 <b>US:</b> 505-1182	

DETACHABLE SIDE RAILS

70.0 cm  
27.5"  
length

Single  
side rail



A stainless steel side rail for standard table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone3 CV	EU: 505-1157
✓ Medstone2 PM	US: 505-1158
✓ Medstone3 PM	
✓ Medstone4 PM	

✓ Included:	! Required:	+ Optional:
One side rail with aluminium double clamps	-	-

Single  
side rail (W)



A stainless steel side rail for wide table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone ERCP	EU: 505-1160
	US: 505-1161

✓ Included:	! Required:	+ Optional:
One side rail with aluminium double clamps	-	-

Pair  
of side rails



Two stainless steel side rails for standard table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone3 CV	EU: 505-1183
✓ Medstone2 PM	US: 505-1184
✓ Medstone3 PM	
✓ Medstone4 PM	

✓ Included:	! Required:	+ Optional:
Two side rails with aluminium double clamps	-	-

Pair  
of side rails (W)



Two stainless steel side rails for wide table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone ERCP	EU: 505-1185
	US: 505-1186

✓ Included:	! Required:	+ Optional:
Two side rails with aluminium double clamps	-	-

Single  
side rail




A stainless steel side rail for standard table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone3 CV	EU: 505-1129
✓ Medstone2 PM	US: 505-1130
✓ Medstone3 PM	
✓ Medstone4 PM	

✓ Included:	! Required:	+ Optional:
One side rail with aluminium double clamps	-	-

Pair  
of side rails

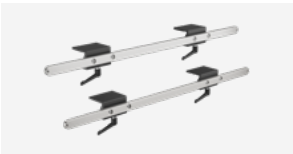

Two stainless steel side rails for standard table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone3 CV	EU: 505-1179
✓ Medstone2 PM	US: 505-1180
✓ Medstone3 PM	
✓ Medstone4 PM	

✓ Included:	! Required:	+ Optional:
Two side rails with aluminium double clamps	-	-

Single  
side rail (W)




A stainless steel side rail for wide table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone ERCP	EU: 505-1155
	US: 505-1156

✓ Included:	! Required:	+ Optional:
One side rail with aluminium double clamps	-	-

Pair  
of side rails (W)

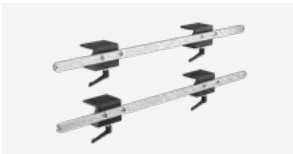

Two stainless steel side rails for wide table tops, to accommodate rail mounted accessories.

Compatibility	Order information
✓ Medstone ERCP	EU: 505-1181
	US: 505-1182

✓ Included:	! Required:	+ Optional:
Two side rails with aluminium double clamps	-	-

80.5 cm  
31.7"  
length



Pair of non symmetrical side rails



Pair of non symmetrical side rails (W)



112.5 cm  
44.3" length

Individual left- and right- stainless steel side rails for accessories on standard table tops.

Individual left- and right- stainless steel side rails for accessories on wide table tops.

Compatibility

Order information

✓

Medstone5 CV

✓

Medstone5 PM

EU:

505-1175

US:

505-1176

Compatibility

Order information

✓

Medstone ERCP

EU:

505-1177

US:

505-1178

Single accessory table clamp



Pair of accessory table clamps



15.2 cm  
6.0" length

One detachable side rail for a light accessory (pan handle, IV-pole) to Medstone3 CV table.

Two detachable side rails used only to attach a catheter tray to Medstone3 CV.

Compatibility

Order information

✓

Medstone3 CV

EU:

505-1119

US:

505-1119

Compatibility

Order information

✓

Medstone3 CV

EU:

505-1121

US\*:

505-1121

\*This pair of clamps is included with a catheter tray for every new Medstone3 CV table (US market only).

- ✓

Included:

Two side rails with aluminium double clamps

!

Required:

-

+

Optional:

-
- ✓

Included:

Two side rails with aluminium double clamps

!

Required:

-

+

Optional:

-

- ✓

Included:

One side rail with aluminium clamp

!

Required:

-

+

Optional:

IV pole  
page 12  
Pan handle  
page 18
- ✓

Included:

Two side rails and two aluminium clamps

!

Required:

-

+

Optional:

Catheter tray  
page 26



Medstone3 ERCP KIT\* (500-1036) or, Medstone3 ERCP (500-1007) with 80.5 cm side rails.

\* Products may not be available in all regions.  
Please contact your STILLE-representative for more information on market availability.

ANESTHESIA ACCESSORIES

Flexible anesthesia screen

Enables hanging cloths or equal material to differentiate the sterile zone from the non-sterile zone. Made in stainless steel and aluminum with polyurethane coating.

- ✔ **Included:**  
One flexible frame with two non-detachable clamps
- ❗ **Required:**  
Fixed or detachable side rails
- ➕ **Optional:**  
-



Compatibility	Order information
✔ Medstone CV	EU: 514-100311
✔ Medstone ERCP	US: 514-100311-US
✔ Medstone PM	

IV pole

Infusion pole with three hooks for side rail fixation. Height is adjustable. Made in stainless steel and aluminum.

- ✔ **Included:**  
One IV pole with a non-detachable clamp
- ❗ **Required:**  
A fixed or detachable side rail
- ➕ **Optional:**  
-



Compatibility	Order information
✔ Medstone CV	EU: 514-100401
✔ Medstone ERCP	US: 514-100401-US
✔ Medstone PM	



Medstone3 PM Rectangular\* (500-1043) with accessory configuration recommended for electropgysiology / pacemaker: 80.5cm side rails (505-1179), X-ray shield loop (535-1772), flexible anesthesia screen (514-100311), IV pole (514-100401), and adjustable arm board (514-50-3SR025).

\* Products may not be available in all regions.  
Please contact your STILLE-representative for more information on market availability.

ARM, LEG & FOOT SUPPORTS

Adjustable arm board

The arm board is designed to support a patient’s arm, ideal for IV. It is very easy to adjust in any desired position with one hand. This arm board is not radiolucent.

- ✔ **Included:**  
One arm board with a non-detachable clamp and an antistatic mattress
- ❗ **Required:**  
Fixed or detachable side rails with 2 clamps (≥ 70.0 cm / 27.5")
- ⊕ **Optional:**  
-



Compatibility	Order information
✔ Medstone CV	EU: 514-50-3SR025
✔ Medstone ERCP	US: 514-50-3SR027
✔ Medstone PM	

Fistula arm board

Radiolucent carbon fiber board with mattress that slides under the table mattress, a fold keeps it from tilting. This accessory provides an unobstructed image.

- ✔ **Included:**  
One carbon fiber board with mattress
- ❗ **Required:**  
-
- ⊕ **Optional:**  
-



Compatibility	Order information
✔ Medstone CV	EU: 505-1118
✔ Medstone PM	US: 505-1118

Dual patient arm board

Radiolucent arm boards, which can rotate horizontally 180 degrees. No rotation lock.

Compatibility	Order information
✔ Medstone5 CV	EU: 505-1110
✔ Medstone PM	US: 505-1110

- ✔ **Included:**  
Two boards in one piece with and mattresses
- ❗ **Required:**  
-
- ⊕ **Optional:**  
-

Combination patient arm board

Radiolucent arm board strongly fixed to the table top and can adjust horizontally 180°.

Compatibility	Order information
✔ Medstone CV	EU: 505-1112
✔ Medstone PM	US: 505-1112

- ✔ **Included:**  
One board with a clamp and mattress
- ❗ **Required:**  
-
- ⊕ **Optional:**  
-

Snap-on arm board

Radiolucent arm board which can be adjusted horizontally 180 degrees.

Compatibility	Order information
✔ Medstone CV	EU: 505-1111
✔ Medstone ERCP	US: Not available
✔ Medstone PM	

- ✔ **Included:**  
One board with a clamp and a mattress
- ❗ **Required:**  
Fixed or detachable side rail
- ⊕ **Optional:**  
-

Basic arm board

Radiolucent arm board with a simple design. Only available in the USA.

Compatibility	Order information
✔ Medstone CV	EU: Not available
✔ Medstone ERCP	US: 512-5000
✔ Medstone PM	


- ✔ **Included:**  
One board, clamp and mattress
- ❗ **Required:**  
Fixed or detachable side rail
- ⊕ **Optional:**  
-





## Arm & leg guard curved

Arm and leg guard intended to support and keep the patient’s arm and leg on the OR-table during surgery and transport. The curve bends outward 5.3 cm (i.e. orthogonal distance between the turning point of the curve and the mattress).

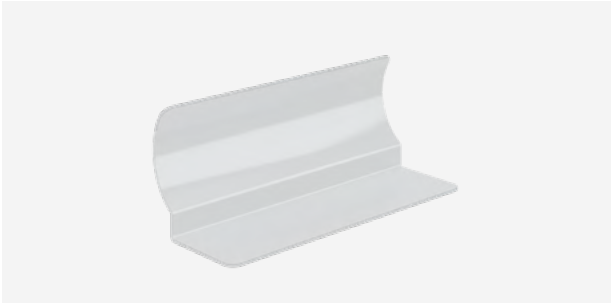- ✔ **Included:**  
One piece of curved guard
- ⓘ **Required:**  
–
- ⊕ **Optional:**  
Cover for Arm & leg guard curved  
EU/US: 514-100553





 Compatibility	 Order information
✔ Medstone CV	EU: 514-100552
✔ Medstone ERCP	US: 514-100552
✔ Medstone PM	

## Arm & leg guard (straight)

Arm and leg guard intended to support and keep the patient’s arm and leg on the OR-table during surgery and transport. The guard slides under the table pad.

- ✔ **Included:**  
One piece of guard
- ⓘ **Required:**  
–
- ⊕ **Optional:**  
Cover for Arm & leg guard straight  
EU/US: 514-100555





 Compatibility	 Order information
✔ Medstone CV	EU: 514-100554
✔ Medstone ERCP	US: 514-100554
✔ Medstone PM	

## Tibia support


The Tibia Support is designed to support the knee or the foot of a patient during knee or leg surgery. Made in stainless steel and includes an antistatic cushion. Cushion diameter is 13.5 cm / 5.3“, length 27 cm / 10.6“.

- ✔ **Included:**  
Metal pole, mattress. The clamp is not included.
- ⓘ **Required:**  
Fixed or detachable side rail  
+ an accessory clamp
- ⊕ **Optional:**  
–




 Compatibility	 Order information
✔ Medstone CV	EU: 514-100367-K
✔ Medstone ERCP	US: Not available
✔ Medstone PM	

## Accessory clamp basic

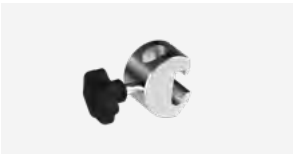

A simple clamp to mount accessories like the tibia support. Made in aluminium.



 Compatibility	 Order information
✔ Medstone CV	EU: 514-100301
✔ Medstone ERCP	US: 514-100301-US
✔ Medstone PM	

- ✔ **Included:**  
One clamp
- ⓘ **Required:**  
–
- ⊕ **Optional:**  
Tibia support (see above)

## Accessory clamp



A robust clamp to mount accessories like the tibia support. Made in stainless steel.

 Compatibility	 Order information
✔ Medstone CV	EU: 514-100305
✔ Medstone ERCP	US: 514-100305-US
✔ Medstone PM	

- ✔ **Included:**  
One clamp
- ⓘ **Required:**  
–
- ⊕ **Optional:**  
Tibia support (see above)

CONTROL DEVICES

Extra pan handle

An additional mushroom switch to move the table top in any horizontal direction with the floating function. Two pan handles can be used simultaneously (one is included with every new Medstone3 CV table).

- ✓ **Included:**  
One pan handle with a non-detachable clamp
- ⓘ **Required:**  
Fixed or detachable side rail
- + **Optional:**  
Accessory table clamp with 15.2 cm / 6" rail – page 11  
EU / US: 505-1119





 Compatibility	 Order information
✓ Medstone3 CV	EU: 505-1171 US: 505-1172

Extension cable

2.5 m / 98.4" extension cable that is connected between the hand control or, a foot control and the control unit. Not for Medstone3 CV. Not available for tables older than 2017.

- ✓ **Included:**  
Cable with connectors
- ⓘ **Required:**  
Hand control (included with every new Medstone table) or, a foot control
- + **Optional:**  
–

 Compatibility	 Order information
✓ Medstone5 CV	EU: 505-1137
✓ Medstone ERCP	US: 505-1137
✓ Medstone PM	

Foot control

Compact multipedal footswitch for unobtrusive placement beneath the table. It has an aluminium body with aluminium or rubber push pedals (depending on which model) and a cable.

- ✓ **Included:**  
One foot switch control with cable and connector
- ⓘ **Required:**  
–
- + **Optional:**  
Extension cable – page 18  
EU / US: 505-1137

Two (2)-movement foot switch control





Maneuvers height adjustment and longitudinal travel.

 Compatibility	 Order information
✓ Medstone2 PM 500-1009 only	EU: 505-1106 US: 505-1106

Foot control for Medstone2 (500-1041, -1011, -1038)



Maneuvers height adjustment and lateral tilt.

 Compatibility	 Order information
✓ Medstone2 ERCP 500-1041	EU: 505-1140 US: 505-1140
✓ Medstone2 PM 500-1011	
✓ Medstone2 PM Rectang. 500-1038	

Foot control for 2 PM (500-1013)



Maneuvers height adjustment and Trendelenburg.

 Compatibility	 Order information
✓ Medstone2 PM 500-1013 only	EU: 505-1141 US: 505-1141

Three (3)-movement  
foot switch control



Maneuvers height adjustment, lateral tilt and Trendelenburg.

Compatibility

- ✓ Medstone3 ERCP 500-1044
- ✓ Medstone3 PM 500-1006
- ✓ Medstone3 PM 500-1039

Order information

EU: 505-1107  
US: 505-1107

Foot control for Medstone3  
(500-1044, -1005, -1043)



Maneuvers height adjustment, longitudinal travel and Trendelenburg.

Compatibility

- ✓ Medstone3 ERCP 500-1044
- ✓ Medstone3 PM 500-1005
- ✓ Medstone3 PM Rectang. 500-1043

Order information

EU: 505-1138  
US: 505-1138

Four (4)-movement  
foot switch control



Maneuvers height adjustment, lateral tilt, Trendelenburg and longitudinal travel.

Compatibility

- ✓ Medstone4 PM 500-1004 only

Order information

EU: 505-1108  
US: 505-1108

Five (5)-movement  
foot switch control



Maneuvers height adjustment, lateral tilt, Trendelenburg, lateral and longitudinal travels.

Compatibility

- ✓ Medstone5 CV 500-1002
- ✓ Medstone5 ERCP 500-1031
- ✓ Medstone5 ERCP KIT 500-1037
- ✓ Medstone5 PM 500-1001

Order information

EU: 505-1109  
US: 505-1109

Foot control for Medstone3  
(500-1007, -1036)



Maneuvers height adjustment, lateral and longitudinal travels.

Compatibility

- ✓ Medstone3 ERCP 500-1007
- ✓ Medstone3 ERCP KIT 500-1036

Order information

EU: 505-1139  
US: 505-1139

A foot control is included with every new Medstone3 CV



Medstone5 PM (500-1001) with X-ray shields kit (535-1770) and foot control (505-1109).  
\* Products may not be available in all regions.  
Please contact your STILLE-representative for more information on market availability.

GENERAL BODY SUPPORTS & BARRIERS

Lateral support

The support can be adjusted in height, longitudinal and lateral, depending on how the patient shall be positioned. Dimensions of the curved cushion 20x20 cm / 7.9" x 7.9".

- ✔ **Included:**  
One metal pole, clamp, curved mattress
- ❗ **Required:**  
Fixed or detachable side rails with 2 clamps (≥ 70.0 cm / 27.5")
- ⊕ **Optional:**  
-



Compatibility	Order information
✔ Medstone CV	EU: 514-50-3SR028
✔ Medstone ERCP	US: 514-50-3SR030
✔ Medstone PM	

Leg / body strap



2.1m / 82.7" strap with velcro fastening, to keep the patient secure to the operating table.

Compatibility	Order information
✔ Medstone CV	EU: 514-100465
✔ Medstone ERCP	US: Not available
✔ Medstone PM	


3 patient straps\* are included with every new Medstone table

\* models may vary according to which region

- ✔ **Included:**  
One patient restraint
- ❗ **Required:**  
-
- ⊕ **Optional:**  
Side rails or side gates - page 22

Pair of side gates

Safety side rails are designed for supporting a patient. They particularly ensure the patient's positioning and comfort during ERCP procedures. They can be easily folded down and when folded to upright position, lock automatically.

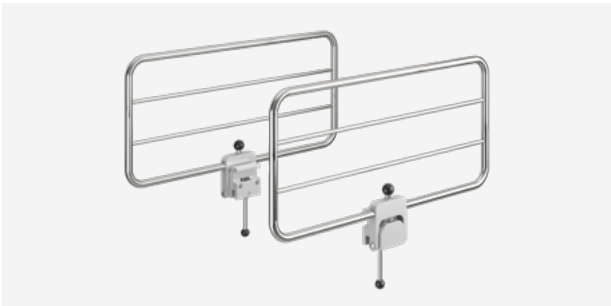

- ✔ **Included:**  
Two side gates, two non-detachable clamps
- ❗ **Required:**  
Fixed or detachable side rails with 2 clamps (≥ 70.0 cm / 27.5")
- ⊕ **Optional:**  
Pair of foldable side gates KIT (two side gates with 2 antistatic covers):  
EU: 514-50-3SR034  
US: 514-50-3SR035  
  
or, single unit of side gates cover:  
EU only: 514-100530-H

Compatibility	Order information
✔ Medstone CV	EU: 514-100530
✔ Medstone ERCP	US: 514-100530-US
✔ Medstone PM	

Patient restraint for standard table top

Padded 3" wide x 70" long for Medstone tables with standard table tops.

Compatibility	Order information
✔ Medstone CV	EU: Not available
✔ Medstone PM	US: 29-034

- ✔ **Included:**  
One padded patient restraint
- ❗ **Required:**  
-
- ⊕ **Optional:**  
Side rails or side gates - page 22

Patient restraint for wide table top

Padded 3" wide x 96" long for Medstone tables with wide table tops.

Compatibility	Order information
✔ Medstone ERCP	EU: Not available
	US: 29-035

- ✔ **Included:**  
One padded patient restraint
- ❗ **Required:**  
-
- ⊕ **Optional:**  
Side rails or side gates - page 22

RADIATION PROTECTION

X-ray protective loop

Radiation shield with lead rubber (0.5 Pb) that protects against scattered radiation. Attach to side rails or directly to the table top with Velcro® strips.



- ✓

**Included:**  
5 protective panels that are fixed to each other
- !

**Required:**  
-
- +

**Optional:**  
Side rail(s)  
X-ray protective loop, grey version  
EU / US: 535-1772  
X-ray protective single shield, grey version  
EU / US: 535-1774  
X-ray protective single shield, blue version  
EU / US: 535-1775



 Compatibility	 Order information
✓ Medstone CV	EU: 535-1773 (blue version)
✓ Medstone ERCP	US: 535-1773 (blue version)
✓ Medstone PM	

X-ray protective kit

Radiation shield with lead rubber (0.5 Pb) that protects against scattered radiation. Attach to the table top with strips of adhesive Velcro® that are placed on the table top. The kit can be used with or without side rails.

- ✓

**Included:**  
2 special panels with cutout for the side rail clamps and 4 straight panels without cutout
- !

**Required:**  
-
- +

**Optional:**  
X-ray protective kit, grey version  
EU / US: 535-1770  
X-ray protective single shield, grey version  
EU / US: 535-1774  
X-ray protective single shield, blue version  
EU / US: 535-1775



 Compatibility	 Order information
✓ Medstone CV	EU: 535-1771 (blue version)
✓ Medstone ERCP	US: 535-1771 (blue version)
✓ Medstone PM	



Medstone5 ERCP\* (500-1031)  
with 80.5 cm side rails (505-1181), X-ray shields (535-1772 and 535-1770) and foot control (505-1109).

*\* Products may not be available in all regions.  
Please contact your STILLE-representative for more information on market availability.*



SURGICAL TABLE EXTENSIONS

Medstone catheter tray

A stainless steel plate that can slide on the side rails of the Medstone3 CV table. The tray / plate dimensions are 60.5 cm / 23.8" width and 65.5 cm / 25.8" length.  
Note: a catheter tray is included with every new Medstone3 CV, in the US only.

- ✓ **Included:**  
One catheter tray and attachment system
- ! **Required:**  
Side rails: fixed rails are included with Medstone3 CV, detachable rails are optional (see below)
- + **Optional:**  
Detachable accessory table clamps with rails to mount the catheter tray at the head end of the table – page 11  
EU / US: 505-1121



Compatibility	Order information
✓ Medstone3 CV	EU: 505-1131 US: 505-1131-US

Adjustable tray

Adjustable tray to place instruments and other devices for the surgical team.

- ✓ **Included:**  
One instrument tray and a clamp
- ! **Required:**  
Fixed or detachable side rail
- + **Optional:**  
-



Compatibility	Order information
✓ Medstone CV	EU: 514-50-3SR036
✓ Medstone ERCP	US: 514-50-3SR037
✓ Medstone PM	

Widening sections

Side extensions can widen the table top to accommodate the larger patient or for procedures involving lateral positioning of the patient, such as for ERCP procedures. Each widening section extends the width of table with 10.4 cm / 4.1".

- ✓ **Included:**  
One or two widening sections (depends on the part number), clamp(s), antistatic mattress(es)
- ! **Required:**  
Fixed or detachable side rail(s) with 2 clamps (≥ 70.0 cm / 27.5")
- + **Optional:**  
300 mm version of this pair of widening sections  
EU: 514-100425  
US: Not available  
455 mm version of this pair of widening sections  
EU: 514-100422  
US: Not available  
500 mm version, single widening section  
EU: 514-100423  
US: Not available

240 mm pair



The shortest widening sections available are: 24 cm / 9.4" length. Two pieces.

Compatibility	Order information
✓ Medstone CV	EU: 514-100421
✓ Medstone ERCP	US: Not available
✓ Medstone PM	

Index

Part number	Accessory (Product name)	Page
29-034	Patient restraint for standard table top	23
29-035	Patient restraint for wide table top	23
505-1106	Two (2) Movement Foot Switch Control	19
505-1107	Three (3) Movement Foot Switch Control	20
505-1108	Four (4) Movement Foot Switch Control	21
505-1109	Five (5) Movement Foot Switch Control	21
505-1110	Dual patient arm board	15
505-1111	Snap-on arm board	15
505-1112	Combination patient arm board	15
505-1118	Fistula arm board	14
505-1119	Single accessory table clamp (Accessory table clamp with 15.2 cm / 6" rail)	11
505-1121	Pair of accessory table clamps (Accessory table clamps – pair)	11
505-1129	Single side rail 80.5 cm / 31.7" length (Medstone detachable side rail EU)	09
505-1130	Single side rail 80.5 cm / 31.7" length (Medstone detachable side rail US)	09
505-1131	Medstone catheter tray	26
505-1131-US	Medstone catheter tray US	26
505-1137	Extension cable (Extension cable for hand control)	18
505-1138	Foot control for Medstone3 (500-1044, -1005, -1043) (Foot Control for 3 PM, 500-1005)	20
505-1139	Foot control for Medstone3 (500-1007, -1036) (Foot Control for 3 PM, 500-1007)	20
505-1140	Foot control for Medstone2 (500-1041, -1011, -1038) (Foot Control for 2 PM, 500-1011)	19
505-1141	Foot Control for 2 PM (500-1013)	19
505-1155	Single side rail (W) 80.5 cm / 31.7" length (Medstone detachable side rail EU Wide)	09
505-1156	Single side rail (W) 80.5 cm / 31.7" length (Medstone detachable side rail US Wide)	09
505-1157	Single side rail 70 cm / 27.5" length (Medstone detachable side rail EU)	08
505-1158	Single side rail 70 cm / 27.5" length (Medstone detachable side rail US)	08
505-1160	Single side rail (W) 70 cm / 27.5" length (Medstone detachable side rail EU Wide)	08
505-1161	Single side rail (W) 70 cm / 27.5" length (Medstone detachable side rail US Wide)	08

Part number	Accessory (Product name)	Page
505-1171	Extra pan handle DIN	18
505-1172	Extra pan handle USA	18
505-1175	Pair of side rails 112.5 cm / 44.3" length (side rails non symmetrical KIT EU)	10
505-1176	Pair of side rails 112.5 cm / 44.3" length (side rails non symmetrical KIT US)	10
505-1177	Pair of side rails (W) 112.5 cm/ 44.3" length (side rails non symmetrical KIT EU Wide)	10
505-1178	Pair of side rails (W) 112.5 cm/ 44.3" length (side rails non symmetrical KIT US Wide)	10
505-1179	Pair of side rails 80.5 cm / 31.7" length (Medstone detachable side rail KIT EU)	09
505-1180	Pair of side rails 80.5 cm / 31.7" length (Medstone detachable side rail KIT US)	09
505-1181	Pair of side rails (W) 80.5 cm / 31.7" length (Medstone detachable side rail EU KIT Wide)	09
505-1182	Pair of side rails (W) 80.5 cm / 31.7" length (Medstone detachable side rail US KIT Wide)	09
505-1183	Pair of side rails 70 cm / 27.5" length (Medstone detachable side rail KIT EU)	08
505-1184	Pair of side rails 70 cm / 27.5" length (Medstone detachable side rail KIT US)	08
505-1185	Pair of side rails (W) 70 cm / 27.5" length (Medstone detachable side rail EU Wide)	08
505-1186	Pair of side rails (W) 70 cm / 27.5" length (Medstone detachable side rail US Wide)	08
505-1188	Vascular package	06
505-1189	ERCP package	07
505-1190	ERCP Package 2	07
512-5000	Basic arm board	15
514-100301	Accessory Clamp basic EU	17
514-100301-US	Accessory Clamp basic US	17
514-100305	Accessory Clamp EU	17
514-100305-US	Accessory Clamp US	17
514-100311	Flexible Anesthesia screen EU	12
514-100311-US	Flexible Anesthesia screen US	12
514-100367-K	Tibia support, short without clamp	17
514-100401	IV pole with clamp EU	12
514-100401-US	IV pole with clamp US	12
514-100421	Widening sections 240mm pair EU	27
514-100422	Widening sections 455mm pair EU	27
514-100423	Widening sections 500mm single EU	27
514-100425	Widening sections 300mm pair EU	27
514-100465	Leg / body strap	23

Part number	Accessory (Product name)	Page
514-100530	Side gates, foldable/ pair EU	22
514-100530-H	Side gates cover single unit	22
514-100530-US	Side gates, foldable/ pair US	22
514-100552	Arm & leg guard curved	16
514-100553	Cover for Arm & leg guard curved	16
514-100554	Arm & leg guard (straight)	16
514-100555	Cover for Arm & leg guard (straight)	16
514-50-3SR025	Adjustable arm board compl. DIN	14
514-50-3SR027	Adjustable arm board compl. US	14
514-50-3SR028	Lateral support complete w clamp EU	22
514-50-3SR030	Lateral support complete w clamp US	22
514-50-3SR034	Side gates, foldable/ pair EU KIT	22
514-50-3SR035	Side gates, foldable/ pair US KIT	22
514-50-3SR036	Adjustable tray with clamp	26
514-50-3SR037	Adjustable tray with clamp US	26
535-1770	X-ray protective kit grey	24
535-1771	X-ray protective kit blue	24
535-1772	Radiation shield base loop grey	24
535-1773	Radiation shield base loop blue	24
535-1774	X-ray protective single shield grey	24
535-1775	X-ray protective single shield blue	24

# STILLE – Surgical perfection. For life.

STILLE is one of the oldest medical device companies in the world. At STILLE, we develop, produce and sell high-quality surgical instruments and surgical imaging tables to physicians worldwide, and have done so for more than 180 years.

STILLE's offering consists of a strong brand and products renowned for their superior quality. We manage our business in two distinct business areas – Surgical Instruments and Surgical Imaging Tables.

## Surgical Instruments

In our Surgical Instruments business area, we have a complete offering of handcrafted, forged instrumentation for all fields of open surgery, with a special emphasis on cardiovascular surgery, and reconstructive and aesthetic plastic surgery. Our instruments feature delicate design, industry

leading durability and a unique feel, which has made STILLE a well-known and highly respected supplier by surgeons throughout the world.

## Surgical Imaging Tables

In our Surgical Imaging Tables business unit, we offer specialty tables for minimally invasive procedures, with pricing, safety and features optimized for both in- and outpatient settings. The STILLE Medstone mobile surgical imaging tables contribute to enhance efficiency in the OR, all while representing valued priced solutions with both quality and reliability.

The STILLE imagiQ2 surgical imaging table with its low-dose capabilities, is optimized for advanced mobile Hybrid OR settings as well as peripheral and interventional procedures.

Learn more at [www.stille.se](http://www.stille.se)

### STILLE AB

Ekbacken 11  
SE-644 30, Torshälla  
Sweden

☎ +46 8 588 580 00

✉ [sales@stille.se](mailto:sales@stille.se)

🌐 [www.stille.se](http://www.stille.se)

### STILLE Surgical Inc.

840 Parkview Blvd  
Lombard, IL 60148  
United States

☎ 1.800.665.1614

✉ [sales.us@stille.se](mailto:sales.us@stille.se)

🌐 [www.stille.se/us](http://www.stille.se/us)



Surgical perfection. For life.



These products may not be available in all regions. Please contact your STILLE-representative for more information on market availability.

This document correctly describe these devices at the time the document was printed. Please always refer to the user manual for any product amendments.

The product images shown are for illustration purposes only and may not be an exact representation of the product.