



Surgical perfection. For life.

High Quality Hospital Grade Surgical Imaging Tables for ERCP with Extra Wide Table Top and 227 kg (500 lbs) Weight Capacity

The STILLE Medstone C-arm tables contribute to enhancing efficiency while representing an affordable solution with quality and reliability provided by Swedish technology.

STILLE Medstone ERCP Surgical Imaging Tables

2, 3- and 5- movement C-arm tables with extra wide table top and optional accessories to secure the patient.

- Rectangular 66 cm (26") wide carbon fiber table top
- Weight capacity 227 kg (500 lbs)
- Motorized movements
- Four hospital grade locking castors
- Metal free imaging 147 cm (58")
- Carbon fiber table top with <1 mm Aluminium Equivalence
- 230Vac, 3,4A, 50 Hz and 120Vac, 6.5A, 60 Hz with 230Vac or 120Vac style line cord
- 5 cm (2") table mattress and 3 patient straps
- Wired hand control

2, 3 or 5 Motorized Movement Options Available:

- Height Adjustment
- Longitudinal travel
- Lateral travel
- Trendelenburg
- Lateral tilt



Photo shows Medstone5 ERCP (500-1031) without side rails

Model	Medstone2 ERCP	Medstone3 ERCP	Medstone5 ERCP
Part number	500-1041*	500-1040*	500-1031*
Table top	66 x 208 cm (26" x 82")	66 x 208 cm (26" x 82")	66 x 208 cm (26" x 82")
Height adjustment	66,3 - 101,6 cm (25.39"-39.37")	67,4 - 102,8 cm (25.39"-40.16")	65.5-101cm (26.78"-39.76")
Lateral tilt	+/- 15°	+/- 15°	+/- 15°
Trendelenburg	x	+/- 15°	+/- 15°
Longitudinal travel	x	x	40 cm (15.75")
Lateral travel	x	x	15 cm (6")

Specification tolerance: $\pm 1^\circ$, $\pm .5"$ (1.27 cm)

*Basic configuration without detachable side rails

Optional Accessories

Part number	Description
505-1189	Medstone ERCP accessory package (includes 505-1181; 514-50-3SR025; 514-100530)
505-1140	Foot control for Medstone2 ERCP
505-1107	Foot control for Medstone3 ERCP
505-1109	Foot control for Medstone5 ERCP
505-1181	Two detachable side rails length 80.5cm (31.69") for Medstone models 500-1040, 500-1041 and 500-1031
514-50-3SR025	Universal adjustable arm board with pad and clamp*
514-100530	Side gates length 61.6 cm (24.25")/ pair with clamps*
514-100311	Anesthesia frame with clamps*
514-100401	IV-pole with clamp*
535-1770	X-ray protective kit grey, 6 single panels; attachment to table top with self-adhesive Velcro © strips
535-1771	X-ray protective kit blue, 6 single panels; attachment to table top with self-adhesive Velcro © strips
535-1772	X-ray protective loop grey, 5 connected panels, attachment to side rails or table top with self-adhesive Velcro © strips*
535-1773	X-ray protective loop blue, 5 connected panels, attachment to side rails or table top with self-adhesive Velcro © strips*


*Accessories attaching to the two detachable side rails 505-1181.



505-1181
Detachable side rails (pair)



514-50-3SR025
Adjustable arm board with pad

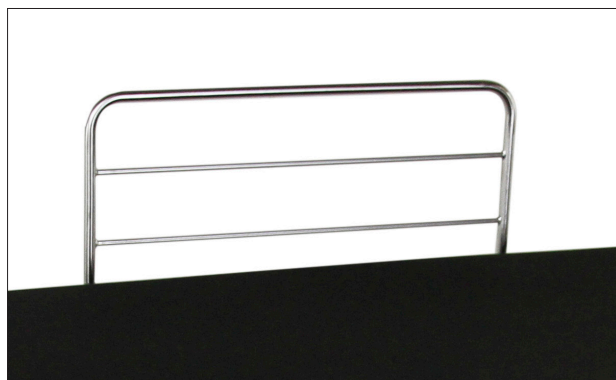

514-100530
Side gates (pair)



535-1772
Radiation shield loop grey



Surgical perfection. For life.

STILLE AB

Ekbacken 11
SE-644 30
Torshälla, Sweden
Phone: +46 8 588 580 00
Fax: +46 8 588 580 05
sales@stille.se

www.stille.se